

# **pUB6/V5-His A, B, and C**

**Catalog no. V250-20**

**Version B**

041202  
25-0251



[www.invitrogen.com](http://www.invitrogen.com)  
[tech\\_service@invitrogen.com](mailto:tech_service@invitrogen.com)

# Table of Contents

Table of Contents.....	iii
General Information.....	iv
<b>Methods.....</b>	<b>1</b>
Overview.....	1
Accessory Products.....	2
Cloning into pUB6/V5-His A, B, and C.....	3
Transfection and Analysis.....	8
Creation of Stable Cell Lines.....	11
<b>Appendix .....</b>	<b>14</b>
Human UbC Promoter .....	14
pUB6/V5-His Vector.....	15
pUB6/V5-His/ <i>lacZ</i> .....	17
Technical Service.....	18
Purchaser Notification .....	20
References.....	21

## General Information

---

### Contents

20 µg each pUB6/V5-His A, B, and C, lyophilized in TE, pH 8.0  
20 µg pUB6/V5-His/lacZ, lyophilized in TE, pH 8.0

---

### Shipping/Storage

Lyophilized plasmids are shipped at room temperature and stored at -20°C.

---

### Product Qualification

pUB6/V5-His A, B, and C and pUB6/V5-His/lacZ are qualified by restriction digest. The table below lists the restriction enzymes and the expected fragments.

Vector	Restriction Enzyme	Expected Fragment
pUB6/V5-His A	<i>Apa</i> I	5462 bp
	<i>Sac</i> II	5462 bp
	<i>Bst</i> E II	No cut
	<i>Stu</i> I	3684 bp, 1778 bp
pUB6/V5-His B	<i>Apa</i> I	5466 bp
	<i>Sac</i> II	4193 bp, 1273 bp
	<i>Bst</i> E II	No cut
	<i>Stu</i> I	3684 bp, 1782 bp
pUB6/V5-His C	<i>Apa</i> I	No cut
	<i>Sac</i> II	5462 bp
	<i>Bst</i> E II	5462 bp
	<i>Stu</i> I	3688 bp, 1774 bp
pUB6/V5-His/lacZ	<i>Apa</i> I	No cut
	<i>Sac</i> II	8509 bp
	<i>Bst</i> E II	8509 bp
	<i>Stu</i> I	4825 bp, 3684 bp

---

# Methods

## Overview

### Introduction

pUB6/V5-His A, B, and C are 5.5 kb vectors designed for overproduction of recombinant proteins in mammalian cell lines. Features of the vectors allow purification and detection of expressed proteins (see pages 15-16 for more information). High-level stable and transient expression can be carried out in most mammalian cells. The vectors contain the following elements:

- Human ubiquitin C promoter (hUbC) for high-level expression across a broad range of species and cell types (Schorpp *et al.*, 1996; Wulff *et al.*, 1990) (see page 14 for more information).
- Three reading frames to facilitate in-frame cloning with a C-terminal peptide encoding the V5 epitope and a polyhistidine (6xHis) metal-binding tag.
- Blasticidin resistance gene (*bsd*) for selection of stable cell lines (see page 11 for more information).
- Episomal replication in cell lines that are latently infected with SV40 or that express the SV40 large T antigen (e.g. COS-1, COS-7).

The control plasmid, pUB6/V5-His/*lacZ*, is included for use as a positive control for transfection, expression, and detection in the cell line of choice.

### Experimental Outline

Use the following outline to clone and express your gene of interest in pUB6/V5-His.

1. Consult the multiple cloning sites described on pages 4-6 to determine which vector (A, B, or C) should be used to clone your gene in frame with the C-terminal V5 epitope and polyhistidine tag.
2. Ligate your insert into the appropriate vector and transform into *E. coli*. Select transformants on 50 to 100 µg/ml ampicillin (or 50 µg/ml blasticidin).
3. Analyze your transformants for the presence of insert by restriction digestion.
4. Select a transformant with the correct restriction pattern and use sequencing to confirm that your gene is cloned in frame with the C-terminal peptide.
5. Transfect your construct into the cell line of choice. Generate a stable cell line, if desired.
6. Test for expression of your recombinant gene by western blot analysis or other functional assay. For antibodies to the V5 epitope or the polyhistidine, C-terminal tag, see the next page.
7. To purify your recombinant protein, you may use metal-chelating resin such as ProBond™. ProBond™ resin is available separately (see next page for ordering information).

## Accessory Products

### Introduction

The products listed below are designed to help you detect and purify your recombinant fusion protein expressed from pUB6/V5-His. In addition, we have a wide variety of mammalian expression vectors, many of which can be utilized with pUB6/V5-His to express and detect multiple proteins in the same cell (see below).

### Antibodies for Detection

If you do not have an antibody to your protein, Invitrogen offers the Anti-V5 antibodies or the Anti-His(C-term) antibodies to detect your recombinant fusion protein. Horseradish peroxidase(HRP)-conjugated antibodies are available for convenient one-step detection.

Antibody	Epitope	Catalog no.
Anti-V5	Detects 14 amino acid epitope derived from the P and V proteins of the paramyxovirus, SV5 (Southern <i>et al.</i> , 1991): GKPIP NPLLGLDST	R960-25
Anti-V5-HRP	same as above	R961-25
Anti-His(C-term)	Detects the C-terminal polyhistidine tag (requires the free carboxyl group for detection) (Lindner <i>et al.</i> , 1997): HHHHHH-COOH	R930-25
Anti-His(C-term)-HRP	same as above	R931-25

### ProBond™ Resin

Ordering information for ProBond™ resin is provided below.

Item	Amount	Catalog no.
ProBond™ Purification System	12 ml precharged ProBond™ resin columns, and buffers for native and denaturing purification	K850-01
ProBond™ Purification Kit with Anti-V5-HRP Antibody	1 Kit	K854-01
ProBond™ Resin	50 ml	R801-01
	150 ml	R801-15

### Expression Vectors

We have a wide variety of mammalian expression vectors utilizing the CMV or EF-1 $\alpha$  promoters. Vectors are available with the Xpress™ (N-terminal), *c-myc* (C-terminal), or V5 (C-terminal) epitopes for detection and either the neomycin, blasticidin, or Zeocin™ resistance genes. All vectors utilize the polyhistidine tag for purification. For more information on the mammalian expression vectors available, see our web site ([www.invitrogen.com](http://www.invitrogen.com)) or call Technical Service (page 18).

# Cloning into pUB6/V5-His A, B, and C

---

## Introduction

Diagrams are provided on pages 4-6 to help you ligate your gene of interest in frame with the coding sequence of the C-terminal peptide. General considerations for transformation are listed below.

---

## General Molecular Biology Techniques

For help with DNA ligations, *E. coli* transformations, restriction enzyme analysis, purification of single-stranded DNA, DNA sequencing, and DNA biochemistry, refer to *Molecular Cloning: A Laboratory Manual* (Sambrook *et al.*, 1989) or *Current Protocols in Molecular Biology* (Ausubel *et al.*, 1994) (See References, page 21).

---

## *E. coli* Strain for Transformation

Many *E. coli* strains are suitable for the propagation of this vector, including TOP10F' (Catalog no. C615-00), JM109, and INVαF' (Catalog no. C658-00). We recommend that you propagate vectors containing inserts in *E. coli* strains that are recombination deficient (*recA*) and endonuclease A deficient (*endA*).

For your convenience, TOP10F' is available as chemically competent or electrocompetent cells from Invitrogen.

Item	Quantity	Catalog no.
Electrocomp™ TOP10F'	5 x 80 µl	C665-55
Ultracomp™ TOP10F' (chemically competent cells)	5 x 300 µl	C665-03
One Shot™ TOP10F' (chemically competent cells)	21 x 50 µl	C3030-03

---

## Transformation Method

You may use any method of your choice for transformation. Chemical transformation is the most convenient for most researchers. Electroporation is the most efficient and the method of choice for large plasmids.

---

## Maintenance of pUB6/V5-His

To propagate and maintain the pUB6/V5-His vectors, resuspend each vector in 20 µl sterile water to prepare a 1 µg/µl stock solution. Store the stock solution at -20°C.

Use this stock solution to transform a *recA*, *endA* *E. coli* strain like TOP10F', DH5α, JM109, or equivalent. Select transformants on LB plates containing 50 to 100 µg/ml ampicillin (or 50 µg/ml blasticidin). Be sure to prepare a glycerol stock of plasmid-containing *E. coli* strain for long-term storage (see page 7).

---

*continued on next page*

## Cloning into pUB6/V5-His A, B, and C, continued

### Cloning Considerations

Your insert should contain a Kozak translation initiation sequence for proper initiation of translation (Kozak, 1987; Kozak, 1991; Kozak, 1990). An example of a Kozak consensus sequence is provided below. Note that other sequences are possible (see references above), but the A at position -3 and the G at position +4 are the most critical (shown in bold). The ATG initiation codon is shown underlined.

ANN**AT**GG

To express your gene as a recombinant fusion protein, you must clone your gene in frame with the C-terminal peptide. The vector is supplied in three reading frames to facilitate cloning. See below and pages 5-6 to develop a cloning strategy.

If you wish to express your protein WITHOUT the C-terminal peptide, be sure to include a stop codon.

### Multiple Cloning Site of pUB6/V5-His A

Below is the multiple cloning site for pUB6/V5-His A. Restriction sites are labeled to indicate the cleavage site. The boxed nucleotides indicate the variable region. Note that there is a stop codon between the *Bam*H I site and the *Bst*X I site. The multiple cloning site has been confirmed by sequencing and functional testing. The complete sequence may be downloaded from our web site ([www.invitrogen.com](http://www.invitrogen.com)) or from Technical Service (see page 18). For more information on the hUbc promoter, see page 14.

```

1121   TTTAGGCACC TTTTGAAATG TAATCATTG GGTCAATATG TAATTTTCAG TGTTAGACTA GTAAATTGTC CGCTAAATTC
                                     UB forward priming site

1201   TGGCCGTTTT TGGCTTTTTT GTTAGAC GAA GCT TGG TAC CGA GCT CGG ATC CAC TAG TCC AGT GTG GTG
                                     Glu Ala Trp Tyr Arg Ala Arg Ile His *** Ser Ser Val Val
                                     Hind III Asp718 I Kpn I Sac I BamH I BstX I*

1270   GAA TTC TGC AGA TAT CCA GCA CAG TGG CGG CCG CTC GAG TCT AGA GGG CCC TTC GAA GGT AAG CCT
                                     EcoR I Pst I EcoR V BstX I* Not I Xho I Xba I Apa I BstB I
                                     Glu Phe Cys Arg Tyr Pro Ala Gln Trp Arg Pro Leu Glu Ser Arg Gly Pro Phe Glu Gly Lys Pro

1336   ATC CCT AAC CCT CTC CTC GGT CTC GAT TCT ACG CGT ACC GGT CAT CAT CAC CAT CAC CAT TGA GT
                                     V5 epitope Age I Polyhistidine region
                                     Ile Pro Asn Pro Leu Leu Gly Leu Asp Ser Thr Arg Thr Gly His His His His His His ***

1401   TTAAACCCGC TGATCAGCCT CGACTGTGCC TTCTAGTTGC CAGCCATCTG TTGTTTGCCC CTCCCCCGTG CCTTCCTTGA
                                     Pme I pcDNA3.1/BGH reverse priming site

1481   CCCTGGAAGG TGCCACTCCC ACTGTCCTTT CCTAATAAAA TGAGGAAATT GCATCGCATT GTCTGAGTAG GTGTCATTCT
                                     BGH polyadenylation signal
  
```

\*Note that there are two *Bst*X I sites in the polylinker.

*continued on next page*

## Cloning into pUB6/V5-His A, B, and C, continued

### Multiple Cloning Site of pUB6/V5-His B

Below is the multiple cloning site for pUB6/V5-His B. Restriction sites are labeled to indicate the cleavage site. The boxed nucleotides indicate the variable region. The multiple cloning site has been confirmed by sequencing and functional testing. The complete sequence may be downloaded from our web site ([www.invitrogen.com](http://www.invitrogen.com)) or from Technical Service (see page 18). For more information on the hUbc promoter, see page 14.

```

1121  TTAGGCACC TTTGAAATG TAATCATTTG GGTCAATATG TAATTTTCAG TGTTAGACTA GTAAATTGTC CGCTAAATTC
                                     UB forward priming site
1201  TGGCCGTTTT TGGCTTTTTT GTTAGACG AAG CTT GGT ACC GAG CTC GGA TCC ACT AGT CCA GTG TGG TGG
                                     Hind III  Asp718 I  Kpn I      Sac I BamH I      BstX I* EcoR I
                                     Lys Leu Gly Thr Glu Leu Gly Ser Thr Ser Pro Val Trp Trp
1271  AAT TCT GCA GAT ATC CAG CAC AGT GGC GGC CGC TCG ACT CTA GAG GGC CCG CGG TTC GAA GGT AAG
                                     Pst I  EcoR V      BstX I*  Not I      Xho I      Xba I      Age I      BstB I
                                     Asn Ser Ala Asp Ile Gln His Ser Gly Gly Arg Ser Ser Leu Glu Gly Pro Arg Phe Glu Gly Lys
1337  CCT ATC CCT AAC CCT CTC CTC GGT CTC GAT TCT ACG CGT ACC GGT CAT CAT CAC CAT CAC CAT TGA
                                     V5 epitope      Age I      Polyhistidine region
                                     Pro Ile Pro Asn Pro Leu Leu Gly Leu Asp Ser Thr Arg Thr Gly His His His His His His ***
1403  GTTTAAAC CCGCTGATCA GCCTCGACTG TGCCTTCTAG TTGCCAGCCA TCTGTTGTTT GCCCCTCCCC CGTGCCTTCC
                                     Pme I      pcDNA3.1/BGH reverse priming site
1481  TTGACCCTGG AAGGTGCCAC TCCCACTGTC CTTTCCTAAT AAAATGAGGA AATTGCATCG CATTGTCTGA GTAGGTGTCA
                                     BGH polyadenylation signal

```

\*Note: that there are two *BstX I* sites in the polylinker.

continued on next page



## Cloning into pUB6/V5-His A, B, and C, continued

### Multiple Cloning Site of pUB6/V5-His C

Below is the multiple cloning site for pUB6/V5-His C. Restriction sites are labeled to indicate the cleavage site. The boxed nucleotides indicate the variable region. The multiple cloning site has been confirmed by sequencing and functional testing. The complete sequence may be downloaded from our web site ([www.invitrogen.com](http://www.invitrogen.com)) or from Technical Service (see page 18). For more information on the hUbC promoter, see page 14.

```

1121   TTTAGGCACC TTTTGAAATG TAATCATTG GGTCAATATG TAATTTTCAG TGTTAGACTA GTAAATTGTC CGCTAAATTC
                                     UB forward priming site
1201   TGGCCGTTTT TGGCTTTTTT GTTAGACGA AGC TTG GTA CCG AGC TCG GAT CCA CTA GTC CAG TGT GGT GGA
                                     Hind III  Asp718 I  Kpn I      Sac I  BamH I                      BstX I*  EcoR I
                                     Ser Leu Val Pro Ser Ser Asp Pro Leu Val Gln Cys Gly Gly
1272   ATT CTG CAG ATA TCC AGC ACA GTG GCG GCC GCT CGA GGT CAC CCA TTC GAA GGT AAG CCT ATC CCT
                                     Pst I  EcoR V      BstX I*  Not I      Xho I      BstE II      BstB I      V5 epitope
                                     Ile Leu Gln Ile Ser Ser Thr Val Ala Ala Ala Arg Gly His Pro Phe Glu Gly Lys Pro Ile Pro
1338   AAC CCT CTC CTC GGT CTC GAT TCT ACG CGT ACC GGT CAT CAT CAC CAT CAC CAT TGA GTTTAA
                                     Age I      Polyhistidine region                      Pme I
                                     Asn Pro Leu Leu Gly Leu Asp Ser Thr Arg Thr Gly His His His His His His ***
                                     pcDNA3.1/BGH Reverse priming site
1401   ACCCGCTGAT CAGCCTCGAC TGTGCCTTCT AGTTGCCAGC CATCTGTTGT TTGCCCTCC CCCGTGCCTT CCTTGACCCT
                                     BGH polyadenylation signal
1481   GGAAGGTGCC ACTCCCACTG TCCTTTCCTA ATAAATGAG GAAATTGCAT CGCATTGTCT GAGTAGGTGT CATTCTATTC

```

Note that there are two *BstX I* sites in the polylinker.

*continued on next page*

## Cloning into pUB6/V5-His A, B, and C, continued

### Transformation of Ligation Mixtures



Transform your ligation mixtures into a competent *recA*, *endA* *E. coli* strain (e.g. TOP10F', DH5 $\alpha$ ) and select on LB plates containing 50-100  $\mu$ g/ml ampicillin or 50  $\mu$ g/ml blasticidin. Select 10-20 clones and analyze for the presence and orientation of your insert.

We recommend that you sequence your construct with the UB Forward and pcDNA3.1/ BGH Reverse primer sequences to confirm that your gene is fused in frame with the V5 epitope and the C-terminal polyhistidine tag. Refer to the diagram on pages 4-6 for the sequence and location of the primer binding sites.

For your convenience, Invitrogen offers a custom primer synthesis service. For more information, refer to our Web site ([www.invitrogen.com](http://www.invitrogen.com)) or contact Technical Service (page 18).

### Preparing a Glycerol Stock

Once you have identified the correct clone, purify the colony and make a glycerol stock for long-term storage. It is also a good idea to keep a DNA stock of your plasmid at -20°C.

- Streak the original colony out for single colonies on an LB plate containing 50  $\mu$ g/ml ampicillin (or 50  $\mu$ g/ml blasticidin). Incubate the plate at 37°C overnight.
- Isolate a single colony and inoculate into 1-2 ml of LB containing 50  $\mu$ g/ml ampicillin (or 50  $\mu$ g/ml blasticidin).
- Grow the culture to mid-log phase ( $OD_{600} = 0.5-0.7$ ).
- Mix 0.85 ml of culture with 0.15 ml of sterile glycerol and transfer to a cryovial.
- Store at -80°C.

# Transfection and Analysis

---

## Introduction

Once you have confirmed that your construct is in the correct orientation and fused to the C-terminal peptide (if desired), then you are ready to transfect your cell line of choice. We recommend that you include the positive control vector and a mock transfection to evaluate your results.

---

## Plasmid Preparation

Plasmid DNA for transfection into eukaryotic cells must be very clean and free from phenol and sodium chloride. Contaminants will kill the cells and salt will interfere with lipids, decreasing transfection efficiency. We recommend isolating plasmid DNA using the S.N.A.P.<sup>™</sup> MiniPrep Kit (10-15 µg DNA, Catalog no. K1900-01), the S.N.A.P.<sup>™</sup> MidiPrep Kit (10-200 µg, Catalog no. K1910-01), or CsCl gradient centrifugation.

---

## Methods of Transfection

For established cell lines (e.g. HeLa), consult original references or the supplier of your cell line for the optimal method of transfection. It is recommended that you follow exactly the protocol for your cell line. Pay particular attention to medium requirements, when to pass the cells, and at what dilution to split the cells. Further information is provided in *Current Protocols in Molecular Biology* (Reference section, page 20).

Methods for transfection include calcium phosphate (Chen and Okayama, 1987; Wigler *et al.*, 1977), lipid-mediated (Felgner *et al.*, 1987; Felgner *et al.*, 1989), and electroporation (Chu *et al.*, 1987; Shigekawa and Dower, 1988). Invitrogen offers the Calcium Phosphate Transfection Kit (K2780-01) and a large variety of transfection reagents for mammalian transfection. For more information on the reagents available, see our Web site ([www.invitrogen.com](http://www.invitrogen.com)) or call Technical Service (see page 18).

---

## Positive Control

pUB6/V5-His/lacZ is provided as a positive control vector for mammalian transfection and expression (see page 17). It may be used to optimize transfection conditions for your cell line. The gene encoding  $\beta$ -galactosidase is expressed in mammalian cells under the control of the human ubiquitin C (hUbC) promoter. A successful transfection will result in  $\beta$ -galactosidase expression that can be easily assayed (see next page).

---

*continued on next page*

## Transfection and Analysis, continued

---

### Assay for $\beta$ -galactosidase Activity

You may assay for  $\beta$ -galactosidase expression by activity assay using cell-free lysates (Miller, 1972) or by staining the cells for activity. Invitrogen offers the  $\beta$ -Gal Assay Kit (Catalog no. K1455-01) and the  $\beta$ -Gal Staining Kit (Catalog no. K1465-01) for fast and easy detection of  $\beta$ -galactosidase expression.

---

### Detection of Fusion Proteins

Several antibodies are available from Invitrogen that can be used to detect expression of your fusion protein from pUB6/V5-His (see page 2).

To detect the fusion protein by Western blot, you will need to prepare a cell lysate from transfected cells. We recommend that you perform a time course to optimize expression of the fusion protein (e.g. 24, 48, 72 hours, etc. after transfection). To lyse cells:

1. Wash cell monolayers ( $\sim 10^6$  cells) once with phosphate-buffered saline (PBS).
  2. Scrape cells into 1 ml PBS and pellet the cells at  $1500 \times g$  for 5 minutes.
  3. Resuspend in 50  $\mu$ l Cell Lysis Buffer (see recipe below). Other lysis buffers are suitable.
  4. Incubate cell suspension at  $37^\circ\text{C}$  for 10 minutes to lyse the cells.
  5. Vortex the cell lysate and centrifuge at  $10,000 \times g$  for 10 minutes to pellet nuclei. Transfer the supernatant to a fresh tube. Assay the supernatant for protein concentration. **Note:** Do not use protein assays utilizing Coomassie Blue or other dyes. NP-40 interferes with the binding of the dye with the protein.
  6. Add SDS-PAGE sample buffer to a final concentration of 1X and boil the sample for 5 minutes.
  7. Load 20  $\mu$ g of lysate onto an SDS-PAGE gel and electrophorese. Use the appropriate percentage of acrylamide to resolve your fusion protein.
- 

### Cell Lysis Buffer

50 mM Tris, pH 7.8  
150 mM NaCl  
1% Nonidet P-40

1. This solution can be prepared from the following common stock solutions. For 100 ml, combine:

1 M Tris base	5 ml
5 M NaCl	3 ml
Nonidet P-40	1 ml
2. Bring the volume up to 90 ml with deionized water and adjust the pH to 7.8 with HCl.
3. Bring the volume up to 100 ml. Store at room temperature.

**Note:** Protease inhibitors may be added at the following concentrations:

1 mM PMSF  
1  $\mu$ g/ml Pepstatin  
1  $\mu$ g/ml Leupeptin

---

*continued on next page*

## Transfection and Analysis, continued



### Note

---

The C-terminal peptide containing the V5 epitope and the polyhistidine region will add approximately 5 kDa to the size of your protein.

---

### Purification

You will need  $5 \times 10^6$  to  $1 \times 10^7$  **transfected** cells for purification of your protein on a 2 ml ProBond™ column (or other metal-chelating column). Refer to the manufacturer's instructions before attempting to purify your fusion protein. To prepare cells for lysis, refer to the protocol on page 13.

---

# Creation of Stable Cell Lines

---

## Introduction

The pUB6/V5-His vectors contain the blasticidin resistance gene (*bsd*) to allow for selection of stable cell lines using blasticidin (Kimura *et al.*, 1994). We recommend that you test the sensitivity of your mammalian host cell to blasticidin as natural resistance varies among cell lines. General information and guidelines are provided below for your convenience.

---


## Blasticidin

Blasticidin S HCl is a nucleoside antibiotic isolated from *Streptomyces griseochromogenes* which inhibits protein synthesis in both prokaryotic and eukaryotic cells (Takeuchi *et al.*, 1958; Yamaguchi *et al.*, 1965). Resistance is conferred by expression of either one of two blasticidin S deaminase genes: *bsd* from *Aspergillus terreus* (Kimura *et al.*, 1994) or *bsr* from *Bacillus cereus* (Izumi *et al.*, 1991). These deaminases convert blasticidin S to a non-toxic deaminohydroxy derivative (Izumi *et al.*, 1991).

---

## Molecular Weight, Formula, and Structure

The formula for blasticidin is  $C_{17}H_{26}N_8O_5 \cdot HCl$ , and the molecular weight is 458.9. The diagram below shows the structure of blasticidin.

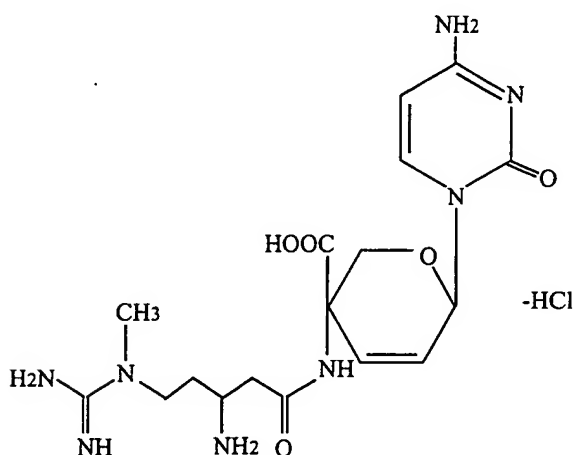

---

## Handling Blasticidin

Always wear gloves, mask, goggles, and protective clothing (e.g. a laboratory coat) when handling blasticidin. Weigh out blasticidin and prepare solutions in a hood.

---

*continued on next page*

## Creation of Stable Cell Lines, continued

### Preparing and Storing Stock Solutions

Blasticidin may be obtained from Invitrogen (Catalog no. R210-01) in 50 mg aliquots. Blasticidin is soluble in water. Water is generally used to prepare stock solutions of 5 to 10 mg/ml.

- Dissolve blasticidin in sterile water and filter-sterilize the solution.
- Aliquot in small volumes suitable for one time use (see last point below) and freeze at -20°C for long-term storage or store at +4°C for short-term storage.
- Aqueous stock solutions are stable for 1-2 weeks at +4°C and 6-8 weeks at -20°C.
- pH of the aqueous solution should not exceed 7 to prevent inactivation of blasticidin.
- Do not subject stock solutions to freeze/thaw cycles (**do not store in a frost-free freezer**).
- Upon thawing, use what you need and discard the unused portion.

### Possible Sites for Linearization

To obtain stable transfectants, you may choose to linearize your vector before transfection. While linearizing your vector may not improve the efficiency of transfection, it increases the chances that the vector does not integrate in a way that disrupts the gene of interest. The table below lists some unique sites that may be used to linearize your construct prior to transformation. **Other restriction sites are possible.** Note that for the enzymes listed below, the cleavage site is indicated for versions A, B, and C of pUB6/V5-His. Be sure that your insert does not contain the restriction enzyme site you wish to use to linearize your vector.

Enzyme	Location	Supplier
<i>Bgl</i> II	Upstream of hUbc promoter	Invitrogen (catalog no. 15213-010)
<i>Bst</i> I 107 I	End of SV40 polyA	AGS*, Fermentas, Takara
<i>Sap</i> I	Backbone	New England Biolabs
<i>Bsp</i> LU I II	Backbone	Boehringer-Mannheim
<i>Alw</i> N I	pMB1 origin	Amersham, New England Biolabs, Life Technologies
<i>Eam</i> I 105 I	Ampicillin gene	AGS*, Fermentas, Takara
<i>Bgl</i> I	Ampicillin gene	Many
<i>Fsp</i> I	Ampicillin gene	Many
<i>Sca</i> I	Ampicillin gene	Invitrogen (catalog no. 15436-017)
<i>Ssp</i> I	Backbone	Invitrogen (catalog no. 15458-011)

*continued on next page*

## Creation of Stable Cell Lines, continued

---

### Selection in Mammalian Cell Lines

To generate a stable cell line expressing your protein, you need to determine the minimum concentration of blasticidin required to kill your untransfected host cell line. Typically, concentrations between 2 and 10 µg/ml blasticidin are sufficient to kill the untransfected host cell line. Test a range of concentrations (see below) to ensure that you determine the minimum concentration necessary for your cell line.

- Seed cells at 20-25% confluency for each time point (~6 time points) and allow the cells to adhere overnight.
  - The next day, substitute culture medium with medium containing varying concentrations of blasticidin (e.g. 0, 1, 3, 5, 7.5, and 10 µg/ml blasticidin).
  - Replenish the selective medium every 3-4 days. Cells sensitive to blasticidin will round up and detach from the plate. Dead cells will accumulate in the medium.
  - Count the number of viable cells at regular intervals to determine the appropriate concentration of blasticidin that prevents growth.
- 

### Selection of Stable Integrants

Once the appropriate concentration of blasticidin is determined, you can generate a stable cell line with your construct.

- Transfect your cells using the optimal protocol for your cell line. Include a sample of untransfected cells as a negative control.
  - 48 hours after transfection, split the cells into fresh medium containing blasticidin at the appropriate concentration for your cell line. Split the cells such that they are no more than 25% confluent.
  - Replenish selective medium every 3-4 days until blasticidin-resistant colonies are detected.
  - Pick and expand colonies.
- 

### Preparation of Cells for Lysis

Use the procedure below to prepare cells for lysis prior to purification of your protein on ProBond™. You will need  $5 \times 10^6$  to  $1 \times 10^7$  cells for purification of your protein on a 2 ml ProBond™ column (see Xpress™ Protein Purification manual).

1. Seed cells (from a stable cell line) in either five T-75 flasks or 2 to 3 T-175 flasks.
  2. Grow the cells in selective medium until they are 80-90% confluent.
  3. Harvest the cells by treating with trypsin-EDTA for 2 to 5 minutes or by scraping the cells in PBS.
  4. Inactivate the trypsin by diluting with fresh medium (if necessary) and transfer the cells to a sterile microcentrifuge tube.
  5. Centrifuge the cells at 1500 x g for 5 minutes. Resuspend the cell pellet in PBS.
  6. Centrifuge the cells at 1500 x g for 5 minutes. You may lyse the cells immediately or freeze in liquid nitrogen and store at -70°C until needed.
- 

### Lysis of Cells

If you are using ProBond™ resin, refer to the ProBond™ Protein Purification manual for details about sample preparation for chromatography.

If you are using other metal-chelating resin, refer to the manufacturer's instruction for recommendations on sample preparation.

---



## Appendix

### Human UbC Promoter

#### Description

The human UbC promoter allows high-level expression of recombinant protein in most mammalian cell lines (Wulff *et al.*, 1990) and in virtually all tissues tested in transgenic mice (Schorpp *et al.*, 1996). The diagram below shows the features of the UbC promoter used in pUB6/V5-His (Nenoi *et al.*, 1996). Features are marked as per Nenoi, *et al.*, 1996.

5' end of UbC promoter

11 GAGATCTGGC CTCCGCGCCG GGTTTTGGCG CCTCCCGCGG GCGCCCCCTT CCTCACGGCG AGCGCTGCCA CGTCAGACGA

Sp 1

91 AGGGCGCAGG AGCGTCCTGA TCCTTCCGCC GGCAGCTCA GGACAGCGGC CCGCTGCTCA TAAGACTCGG CCTTAGAAC

171 CCAGTATCAG CAGAAGGACA TTTTAGGACG GGAATTGGGT GACTCTAGGG CACTGGTTTT CTTCCAGAG AGCGGAACAG

Sp 1 TATA box

251 GCGAGGAAAA GTAGTCCCTT CTCGGCGATT CTGCGGAGGG ATCTCCGTG GCGCGTGAAC GCCGATGATT ATATAAGGAC

Start of Transcription

331 GCGCCGGGTG TGGCACAGCT AGTTCCTGCG CAGCCGGGAT TTGGGTCGCG GTTCTTGTTT GTGGATCGCT GTGATCGTCA

Exon 1

5' end of Intron 1

411 CTTGGTGAGT AGCGGGCTGC TGGGCTGGCC GGGGCTTTTC TGGCCGCCG GCCGCTCGGT GGGACGGAAG CGTGTGGAGA

491 GACCGCCAAG GGCTGTAGTC TGGGTCCGCG AGCAAGGTTG CCCTGAACTG GGGGTTGGGG GGAGCGCAGC AAAATGGCGG

571 CTGTTCCCGA GTCTTGAATG GAAGACGCTT GTGAGGCGGG CTGTGAGGTC GTTGAAACAA GGTGGGGGGC ATGGTGGGCG

651 GCAAGAACC AAGGTCTTGA GGCCTTCGCT AATCGGGAA AGCTCTTATT CGGGTGAGAT GGGCTGGGGC ACCATCTGGG

731 GACCTGACG TGAAGTTTGT CACTGACTGG AGAACTCGGT TTGTCGTCTG TTGCGGGGGC GGCAGTTATG CGGTGCCGTT

811 GGGCAGTGCA CCCGTACCTT TGGGAGCGCG CGCCCTCGTC GTGTCGTGAC GTCACCCGTT CTGTTGGCTT ATAATGCAGG

891 GTGGGGCCAC CTGCCGGTAG GTGTGCGGTA GGCTTTTCTC CGTCGCAGGA CGCAGGGTTC GGGCCTAGGG TAGGCTCTCC

971 TGAATCGACA GCGCGCGGAC CTCTGGTGAG GGGAGGGATA AGTGAGGCGT CAGTTTCTTT GGTGCGTTTT ATGTACCTAT

1051 CTTCTTAAGT AGCTGAAGCT CCGGTTTTGA ACTATGCGCT CGGGGTTGGC GAGTGTGTTT TGTGAAGTTT TTTAGGCACC

UB Forward priming site

1131 TTTTGAAATG TAATCATTTG GGTCAATATG TAATTTTCAG TGTTAGACTA GTAAATTGTC CGCTAAATTC TGGCCGTTTT

3' end of Intron 1


1211 TGGCTTTTTT GTTAGACGAA GCTTGG....

5' end of Exon 2

## pUB6/V5-His Vector

### Map of pUB6/ V5-His

The figure below summarizes the features of the pUB6/V5-His vectors. The sequences for pUB6/V5-His A, B, and C are available for downloading from our Web site (<http://www.invitrogen.com>) or from Technical Service (see page 18).

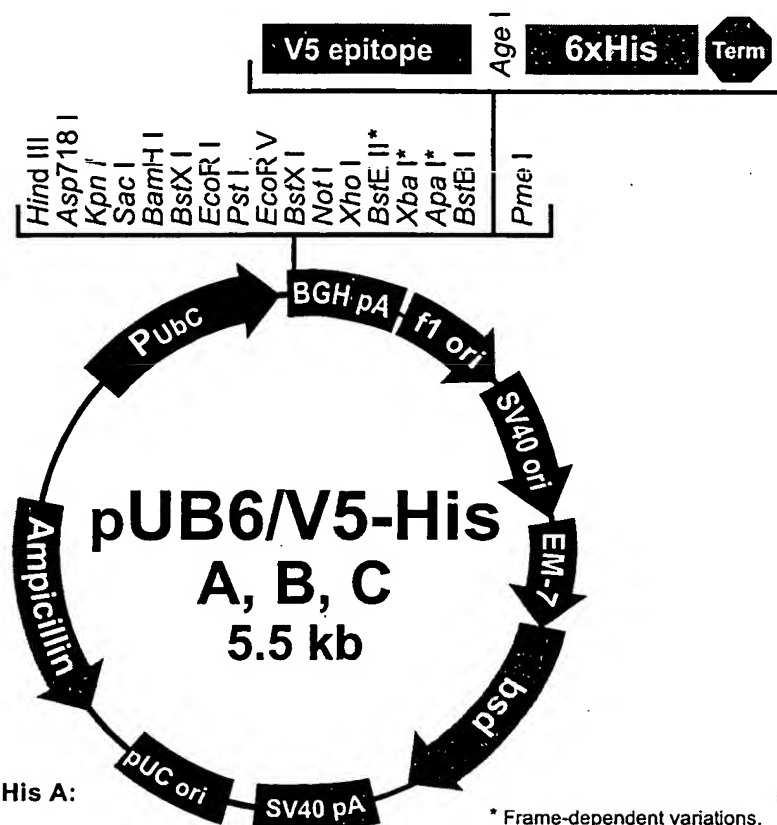

**Comments for pUB6/V5-His A:**  
5463 nucleotides

UbC promoter: bases 18-1227  
 UB forward priming site: bases 1167-1188  
 Multiple cloning site: bases 1229-1326  
 V5 epitope: bases 1327-1368  
 Polyhistidine tag: bases 1378-1395  
 pcDNA3.1/BGH reverse priming site: bases 1418-1435  
 BGH polyadenylation signal: bases 1421-1648  
 f1 origin: bases 1694-2122  
 SV40 promoter and origin: bases 2149-2458  
 EM-7 promoter: bases 2506-2561  
 Blastidicin resistance gene: bases 2580-2978  
 SV40 polyadenylation signal: bases 3136-3266  
 pUC origin: bases 3649-4322  
 Ampicillin resistance gene: bases 4467-5327

\* Frame-dependent variations.  
BstE II is only found in version C.  
In addition, there are no Xba I  
or Apa I sites in version C.

*continued on next page*

## pUB6/V5-His Vector, continued

### Features of pUB6/V5-His

pUB6/V5-His A (5463 bp), pUB6/V5-His B (5467 bp), and pUB6/V5-His C (5459 bp) contain the following elements. All features have been functionally tested.

Feature	Benefit
Human ubiquitin C (hUbC) promoter	Permits overexpression of your recombinant protein in a broad range of mammalian cell types (Hershko and Ciechanover, 1982; Wulff <i>et al.</i> , 1990; Schorpp <i>et al.</i> , 1996)
Multiple cloning site in three reading frames	Allows insertion of your gene and facilitates cloning in frame with the V5 epitope and the C-terminal polyhistidine tag
V5 epitope (Gly-Lys-Pro-Ile-Pro-Asn-Pro-Leu-Leu-Gly-Leu-Asp-Ser-Thr)	Allows detection of your recombinant protein with the Anti-V5 Antibody (Catalog no. R960-25) or Anti-V5-HRP Antibody (Catalog no. R961-25) (Southern <i>et al.</i> , 1991)
C-terminal polyhistidine tag	Permits purification of your recombinant protein on metal-chelating resin such as ProBond™ (R801-01)  In addition, the C-terminal polyhistidine tag is the epitope for the Anti-His(C-term) Antibody (Catalog no. R930-25) and the Anti-His(C-term)-HRP Antibody (Catalog no. R931-25) (Lindner <i>et al.</i> , 1997)
Bovine growth hormone (BGH) polyadenylation signal	Efficient transcription termination and polyadenylation of mRNA (Goodwin and Rottman, 1992)
f1 origin	Allows rescue of single-stranded DNA
SV40 early promoter and origin	Allows efficient, high-level expression of the blasticidin resistance gene and episomal replication in cells expressing the SV40 large T antigen
EM-7 promoter	For expression of the blasticidin resistance gene in <i>E. coli</i> .
Blasticidin resistance gene ( <i>bsd</i> )	Selection of stable transfectants in mammalian cells (Kimura <i>et al.</i> , 1994)
SV40 polyadenylation signal	Efficient transcription termination and polyadenylation of mRNA
pUC origin	High-copy number replication and growth in <i>E. coli</i>
Ampicillin resistance gene ( $\beta$ -lactamase)	Selection of vector in <i>E. coli</i>


## pUB6/V5-His/lacZ

### Description

pUB6/V5-His/lacZ is a 8510 bp control vector containing the gene for  $\beta$ -galactosidase. This vector was constructed by ligating a 3190 bp *Hind* III-*Age* I fragment containing the *lacZ* gene and the V5 epitope to a 5320 bp *Hind* III-*Age* I fragment containing the hUbC promoter, polyhistidine tag and blasticidin resistance gene from pUB6/V5-His A.

### Map of Control Vector

The figure below summarizes the features of the pUB6/V5-His/lacZ vector. The complete nucleotide sequence for pUB6/V5-His/lacZ is available for downloading from our Web site (<http://www.invitrogen.com>) or by contacting Technical Service. (see next page).

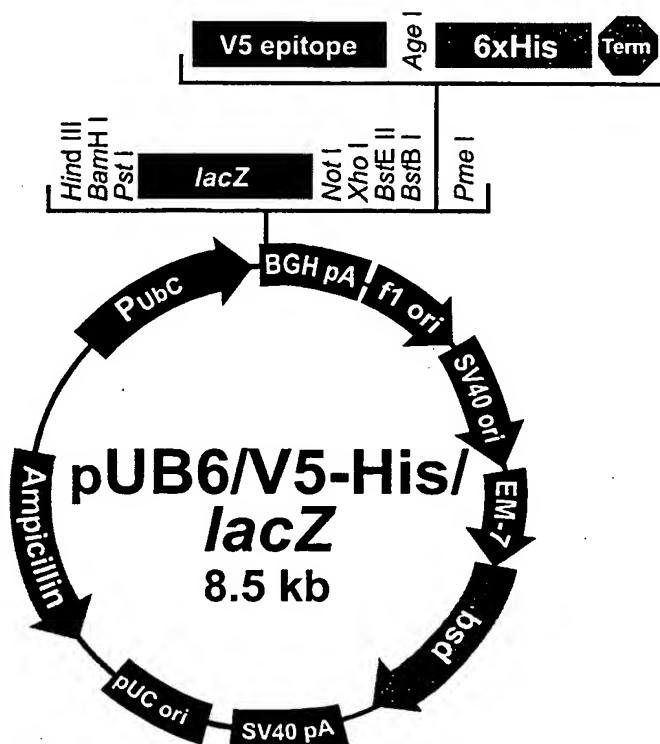

#### Comments for pUB6/V5-His/lacZ: 8510 nucleotides

UbC promoter: bases 18-1227  
UB forward priming site: bases 1167-1188  
*LacZ* ORF: bases 1290-4346  
V5 epitope: bases 4374-4415  
Polyhistidine tag: bases 4425-4442  
pcDNA3.1/BGH reverse priming site: bases 4465-4482  
BGH polyadenylation signal: bases 4468-4695  
f1 origin: bases 4741-5169  
SV40 promoter and origin: bases 5196-5505  
EM-7 promoter: bases 5553-5608  
Blasticidin resistance gene: bases 5627-6025  
SV40 polyadenylation signal: bases 6183-6313  
pUC origin: bases 6696-7369  
Ampicillin resistance gene: bases 7514-8374

## Technical Service

### World Wide Web



Visit the Invitrogen Web Resource using your World Wide Web browser. At the site, you can:

- Get the scoop on our hot new products and special product offers
- View and download vector maps and sequences
- Download manuals in Adobe® Acrobat® (PDF) format
- Explore our catalog with full color graphics
- Obtain citations for Invitrogen products
- Request catalog and product literature

Once connected to the Internet, launch your Web browser (Internet Explorer 5.0 or newer or Netscape 4.0 or newer), then enter the following location (or URL):

**<http://www.invitrogen.com>**

...and the program will connect directly. Click on underlined text or outlined graphics to explore. Don't forget to put a bookmark at our site for easy reference!

### Contact Us

For more information or technical assistance, call, write, fax, or email. Additional international offices are listed on our Web page ([www.invitrogen.com](http://www.invitrogen.com)).

#### Corporate Headquarters:

Invitrogen Corporation  
1600 Faraday Avenue  
Carlsbad, CA 92008  
USA

Tel: 1 760 603 7200

Tel (Toll Free): 1 800 955 6288

Fax: 1 760 602 6500

E-mail:

[tech\\_service@invitrogen.com](mailto:tech_service@invitrogen.com)

#### Japanese Headquarters:

Invitrogen Japan K.K.  
Nihonbashi Hama-Cho Park  
Bldg. 4F  
2-35-4, Hama-Cho, Nihonbashi

Tel: 81 3 3663 7972

Fax: 81 3 3663 8242

E-mail: [jpinfo@invitrogen.com](mailto:jpinfo@invitrogen.com)

#### European Headquarters:

Invitrogen Ltd  
3 Fountain Drive  
Inchinnan Business Park  
Paisley PA4 9RF, UK

Tel: +44 (0) 141 814 6100

Tel (Toll Free): 0800 5345 5345

Fax: +44 (0) 141 814 6287

E-mail: [eurotech@invitrogen.com](mailto:eurotech@invitrogen.com)

*continued on next page*

## Technical Service, continued

---

### Limited Warranty

Invitrogen is committed to providing our customers with high-quality goods and services. Our goal is to ensure that every customer is 100% satisfied with our products and our service. If you should have any questions or concerns about an Invitrogen product or service, please contact our Technical Service Representatives.

Invitrogen warrants that all of its products will perform according to the specifications stated on the certificate of analysis. The company will replace, free of charge, any product that does not meet those specifications. This warranty limits Invitrogen Corporation's liability only to the cost of the product. No warranty is granted for products beyond their listed expiration date. No warranty is applicable unless all product components are stored in accordance with instructions. Invitrogen reserves the right to select the method(s) used to analyze a product unless Invitrogen agrees to a specified method in writing prior to acceptance of the order.

Invitrogen makes every effort to ensure the accuracy of its publications, but realizes that the occasional typographical or other error is inevitable. Therefore Invitrogen makes no warranty of any kind regarding the contents of any publications or documentation. If you discover an error in any of our publications, please report it to our Technical Service Representatives.

Invitrogen assumes no responsibility or liability for any special, incidental, indirect or consequential loss or damage whatsoever. The above limited warranty is sole and exclusive. No other warranty is made, whether expressed or implied, including any warranty of merchantability or fitness for a particular purpose.

---

## Purchaser Notification

<b>Introduction</b>	Use of the pUB6/V5-His A, B, and C is covered under the licenses detailed below.
<b>Limited Use Label License No. 22: Vectors and Clones containing sequences coding for Histidine Hexamer</b>	<p>This product is licensed under U.S. Patent Nos. 5,284,933 and 5,310,663 and foreign equivalents from Hoffmann-LaRoche, Inc., Nutley, NJ and/or Hoffmann-LaRoche Ltd., Basel, Switzerland and is provided only for use in research.</p> <p>Information about licenses for commercial use is available from:</p> <p>QIAGEN GmbH Max-Volmer-Str. 4 D-40724 Hilden Germany</p>
<b>Limited Use Label License No. 49: BGH Polyadenylation Signal</b>	<p>The bovine growth hormone (BGH) polyadenylation sequence is licensed under U.S. Patent No. 5,122,458 for research purposes only. "Research purposes" means uses directed to the identification of useful recombinant proteins and the investigation of the recombinant expression of proteins, which uses shall in no event include any of the following: Any use in humans of a CLAIMED DNA or CLAIMED CELL; Any use in human of protein or other substance expressed or made at any stage of its production with the use of a CLAIMED DNA or a CLAIMED CELL; Any use in which a CLAIMED DNA or CLAIMED CELL would be sold or transferred to another party other than Invitrogen, its affiliate, or its sublicensee; Any use in connection with the expression or production of a product intended for sale or commercial use; or Any use for drug screening or drug development.</p> <p>Inquiries for commercial use should be directed to: Bennett Cohen, Ph.D., Research Corporation Technologies, 101 North Wilmot Road, Suite 600, Tucson, AZ 85711-3335, Tel: 1-520-748-4400 Fax: 1-520-748-0025.</p>
<b>Limited Use Label License No. 51: Blasticidin and the Blasticidin Selection Marker</b>	<p>Blasticidin and the blasticidin resistance gene (<i>bsd</i>) are sold under patent license and may be used for <b>research purposes only</b>. Inquiries for commercial use should be directed to: Kaken Pharmaceutical Company, Ltd., Bunkyo Green Court, Center Office Building, 19-20 Fl, 28-8 Honkomagome 2-chome, Bunkyo-ku, Tokyo 113-8650, Japan.</p> <p>Tel: 81 3-5977-5008; Fax: 81 3-5977-5008.</p>
<b>Product User Registration Card</b>	<p>Please complete and return the enclosed Product User Registration Card for pUB6/V5-His A, B, and C that you purchase. This will serve as a record of your purchase and will allow Invitrogen to provide you with critical product updates. The agreement outlined above becomes effective upon our receipt of your Product User Registration Card or 10 days following the sale of pUB6/V5-His A, B, and C to you. Use of the vectors at any time results in immediate obligation to the terms and conditions stated in this Agreement.</p>

## References

---

- Ausubel, F. M., Brent, R., Kingston, R. E., Moore, D. D., Seidman, J. G., Smith, J. A., and Struhl, K. (1994). *Current Protocols in Molecular Biology* (New York: Greene Publishing Associates and Wiley-Interscience).
- Chen, C., and Okayama, H. (1987). High-Efficiency Transformation of Mammalian Cells by Plasmid DNA. *Mol. Cell. Biol.* **7**, 2745-2752.
- Chu, G., Hayakawa, H., and Berg, P. (1987). Electroporation for the Efficient Transfection of Mammalian Cells with DNA. *Nuc. Acids Res.* **15**, 1311-1326.
- Felgner, P. L., Gadek, T. R., Holm, M., Roman, R., Chan, H. W., Wenz, M., Northrop, J. P., Ringold, G. M., and Danielsen, M. (1987). Lipofectin: A Highly Efficient, Lipid-mediated DNA-transfection Procedure. *Proc. Natl. Acad. Sci. USA* **84**, 7413-7417.
- Felgner, P. L., Holm, M., and Chan, H. (1989). Cationic Liposome Mediated Transfection. *Proc. West. Pharmacol. Soc.* **32**, 115-121.
- Goodwin, E. C., and Rottman, F. M. (1992). The 3'-Flanking Sequence of the Bovine Growth Hormone Gene Contains Novel Elements Required for Efficient and Accurate Polyadenylation. *J. Biol. Chem.* **267**, 16330-16334.
- Hershko, A., and Ciechanover, A. (1982). Mechanisms of Intracellular Protein Breakdown. *Ann. Rev. Biochem.* **51**, 335-364.
- Izumi, M., Miyazawa, H., Kamakura, T., Yamaguchi, I., Endo, T., and Hanaoka, F. (1991). Blasticidin S-Resistance Gene (*bsr*): A Novel Selectable Marker for Mammalian Cells. *Exp. Cell Res.* **197**, 229-233.
- Kimura, M., Takatsuki, A., and Yamaguchi, I. (1994). Blasticidin S Deaminase Gene from *Aspergillus terreus* (*BSD*): A New Drug Resistance Gene for Transfection of Mammalian Cells. *Biochim. Biophys. Acta* **1219**, 653-659.
- Kozak, M. (1987). An Analysis of 5'-Noncoding Sequences from 699 Vertebrate Messenger RNAs. *Nuc. Acids Res.* **15**, 8125-8148.
- Kozak, M. (1991). An Analysis of Vertebrate mRNA Sequences: Intimations of Translational Control. *J. Cell Biol.* **115**, 887-903.
- Kozak, M. (1990). Downstream Secondary Structure Facilitates Recognition of Initiator Codons by Eukaryotic Ribosomes. *Proc. Natl. Acad. Sci. USA* **87**, 8301-8305.
- Lindner, P., Bauer, K., Krebber, A., Nieba, L., Kremmer, E., Krebber, C., Honegger, A., Klinger, B., Mocikat, R., and Pluckthun, A. (1997). Specific Detection of His-tagged Proteins With Recombinant Anti-His Tag scFv-Phosphatase or scFv-Phage Fusions. *BioTechniques* **22**, 140-149.
- Miller, J. H. (1972). *Experiments in Molecular Genetics* (Cold Spring Harbor, New York: Cold Spring Harbor Laboratory).
- Nenoi, M., Mita, K., Ichimura, S., Cartwright, I. L., Takahashi, E., Yamaguchi, M., and Tsuji, H. (1996). Heterogeneous Structure of the Polyubiquitin Gene UbC of HeLa S3 Cells. *Gene* **175**, 179-185.
- Sambrook, J., Fritsch, E. F., and Maniatis, T. (1989). *Molecular Cloning: A Laboratory Manual*, Second Edition (Plainview, New York: Cold Spring Harbor Laboratory Press).

---

*continued on next page*



## References, continued

---

- Schorpp, M., Jäger, R., Schellander, K., Schenkel, J., Wagner, E. F., Weiher, H., and Angel, P. (1996). The Human Ubiquitin C Promoter Directs High Ubiquitous Expression of Transgenes in Mice. *Nuc. Acids Res.* 24, 1787-1788.
- Shigekawa, K., and Dower, W. J. (1988). Electroporation of Eukaryotes and Prokaryotes: A General Approach to the Introduction of Macromolecules into Cells. *BioTechniques* 6, 742-751.
- Southern, J. A., Young, D. F., Heaney, F., Baumgartner, W., and Randall, R. E. (1991). Identification of an Epitope on the P and V Proteins of Simian Virus 5 That Distinguishes Between Two Isolates with Different Biological Characteristics. *J. Gen. Virol.* 72, 1551-1557.
- Takeuchi, S., Hirayama, K., Ueda, K., Sakai, H., and Yonehara, H. (1958). Blastcidin S, A New Antibiotic. *The Journal of Antibiotics, Series A* 11, 1-5.
- Wigler, M., Silverstein, S., Lee, L.-S., Pellicer, A., Cheng, Y.-C., and Axel, R. (1977). Transfer of Purified Herpes Virus Thymidine Kinase Gene to Cultured Mouse Cells. *Cell* 11, 223-232.
- Wulff, B. S., O'Hare, M. M., Boel, E., Theill, L. E., and Schwartz, T. W. (1990). Partial Processing of the Neuropeptide Y Precursor in Transfected CHO Cells. *FEBS Lett.* 261, 101-105.
- Yamaguchi, H., Yamamoto, C., and Tanaka, N. (1965). Inhibition of Protein Synthesis by Blastcidin S. I. Studies with Cell-free Systems from Bacterial and Mammalian Cells. *J. Biochem. (Tokyo)* 57, 667-677.

---

©1998-2002 Invitrogen Corporation. All rights reserved.